



Mortuary Refrigeration



Small Door End Opening Modular Systems
Intelligently Designed From The Inside Out

ANATOMICAL PATHOLOGY | MEDICAL EXAMINER/MORGUE | ANATOMY LABS | NECROPSY

Mopec Delivers More

SMALL DOOR END OPENING MODULAR SYSTEM FEATURES

- Quality construction - state-of-the-art manufacturing facility
- High density insulation - Foamed in place polyurethane insulation
- Temperature control - 20°C (-4F) for freezer, 4°C (39.3F) for refrigeration
- Installation available - Factory trained installation team

AVAILABLE FOR QUICK SHIP

Mopec morgue refrigeration units can be sized to suit your facility needs, no matter the size. Small door end opening modular systems are compatible with either telescoping or conveyor tray storage and are available for immediate shipment.

VERSATILE CONFIGURATION



KM201



KM202



KM203



KM301



KM302




KM303



SUPERIOR DOOR DESIGN

For guaranteed strength, Mopec provides quality flush doors complete with patented extruded anodized aluminum frames. Each unit's interior and exterior is clad of stainless steel that makes cleaning virtually effortless. The magnetic rubber gaskets give the tightest seal, while a heated perimeter and strike surface eliminate condensation.

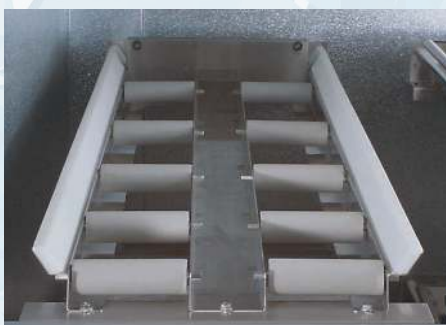

CONVEYOR TRAY STORAGE

The conveyor tray system is permanently located within the cooler and fully supports the body tray. When the door is open, accessibility is limited to use of a lift with the same conveyor system. The tray can then be rolled out of the cooler and onto the conveyor lift system. The conveyor mechanism system has a split wheel design which will accommodate body trays with drain holes.



TELESCOPING TRAY STORAGE


A telescoping mechanism fully supports the body tray. When the door is open and the end of the telescoping mechanism is fully extended, the tray is in full head-to-toe view. The cadaver can be removed alone or with the body tray.

Sanisteel®

Antibacterial Pre-Coated Steel



Sanisteel® is a new film coating with antibacterial properties intended to become the point of reference for the pre-coated steel sector. Sanisteel® activates its antibacterial characteristics in the presence of environmental conditions which favor the proliferation of bacteria, fungi and algae. It has a dual action, inhibits the growth of germs and eliminates those who might be present. It is a focused solution to maintain a long lasting protection and hygiene even in those situations where normal cleaning may not be guaranteed.

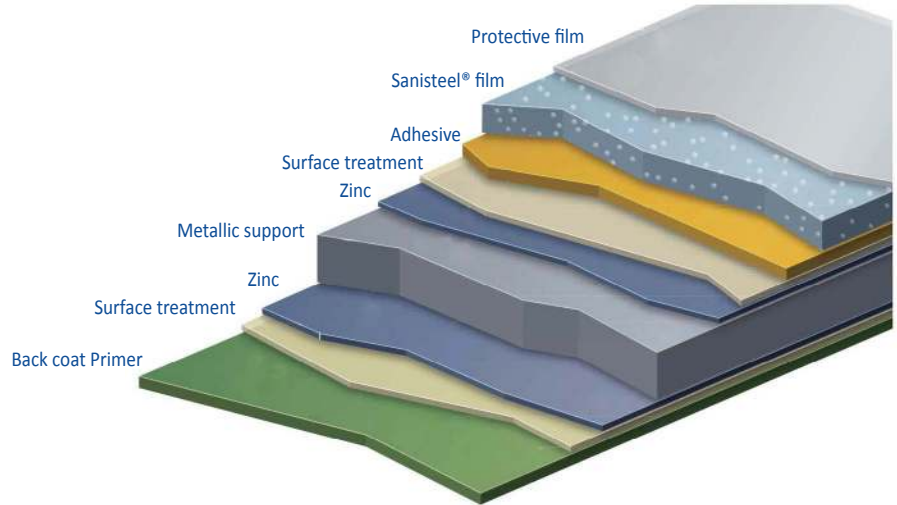

COMPATIBLE ACCESSORIES WITH MOPEC REFRIGERATION SYSTEMS



TRANSPORT We offer body tray carriers, cadaver carriers and cadaver covers.



TRAYS We offer body trays in pediatric, normal, large and bariatric sizes.



BODY BAGS We offer a wide assortment of body bags including pediatric, bariatric, water recovery, infection containment and cremation safe



LIFTS We offer conveyor, pallet, low profile, hydraulic and overhead lifts.

REQUEST A CONSULATION

For product information and design assistance, call: 800.362.8491 or Visit us at www.mopec.com

About Mopec



Mopec is committed to designing and manufacturing equipment and accessories to enhance your facility's safety, flexibility and productivity. Our technical prowess stems from a commitment to quality craftsmanship and high-touch service. Mopec is recognized globally for giving clients an unrivaled experience including design, manufacturing, installation and post-sales support. The industries we serve include anatomic pathology, medical examiner/morgue, anatomy labs, and necropsy.



INDUSTRIES WE SERVE



Anatomic Pathology



Medical Examiner/Morgue



Anatomy Labs



Necropsy

VISIT WWW.MOPEC.COM FOR MORE INFORMATION

21750 Coolidge Highway • Oak Park, MI 48237 USA
1-800-362-8491 • +1 248-291-2040 • Email: info@mopec.com

Equipment shown may include optional features. Please visit our website for more information.

Copyright ©Mopec 2016