

# Mopec *Maestro*

GROSSING STATION



EQUIP YOUR LABORATORY WITH  
**THE SAFEST & MOST PERSONALIZED GROSSING STATION**

# Mopec Maestro<sup>®</sup> Grossing Station

*Mopec understands that each pathology professional has their own workspace preferences, and each laboratory holds their own safety, ergonomic and workflow requirements. With the Maestro's distinctive Ensemble Pegboard™, dozens of accessories and flexible safety options, every station can be quickly personalized to suit each individual.*

## PREMIUM WORKSURFACE

### Cutting & rinsing surface

Each station features two removable panels with one perforated plate and one specimen cutting board with removable legs. The panels are capable of sliding across set rails for workspace flexibility and easy cleaning. The sink can also be adjusted for left or right handed users.

### Height adjustability

Maestro workstations are ADA compliant, and they adjust to accommodate professionals of varying heights.



## INTUITIVE CONTROLS & DEVICE MANAGEMENT

### Easy-to-use touchscreen

With a quick tap, all Maestro features can be controlled using the new, integrated touchscreen.

### Wire containment

Four outlets are provided on each bench with two USB compatible and two GFCI. Cable covers are incorporated into the peg board and on the side of each station, therefore eliminating the look of messy cords.

## TOP-NOTCH VENTILATION & SAFETY

### Front air system

This patented, on-demand technology pushes air from the front of the cutting surface through the back of the area, creating a laminar flow pattern. This continuous ventilation protects the user from dangerous formalin exposure.

### Safety notifications

On-screen warnings alert users of possible harmful conditions involving air flow and optional formalin collection and dispensing systems.



# Harmonize Your Workspace with Optional Accessories

## ENSEMBLE PEGBOARD™ ACCESSORIES

Movable options for precise workstation personalization. Add-on pegboard wings allow for more storage space.

- Organizer Bins
- Requisition File Holder
- Bar Code Holster
- Shelves
- Magnetic Tool Holder
- Staining Dye Holder
- Specimen Scale Shelf
- Towel & Glove Box Holders
- SaniPath™ Wipe Dispenser
- Safety Glasses Holder
- Face Shield Holster
- Monitor Mount

MULTIPLE MONITOR MOUNTING OPTIONS FOR IMPROVED ERGONOMICS AND REDUCED NECK STRAIN.



## OPTIONAL DRAWER SYSTEM

Storage options for quick reach accessibility



- Side Drawer
- Middle Drawer
- Small Over-the-Sink Drawer
- Instrument Soaking Inserts
- Keyboard Tray Insert
- Tissue Marking Dyes Insert



## OTHER OPTIONS & ACCESSORIES

Additions for an enhanced workspace

- Cassette Labeling Shelf
- Side Monitor Arm
- Spatter-Free Formalin Dispenser
- Ventilated Trash
- CPU Bracket
- PathCam™ Imaging System
- Foot Pedal Faucet Control
- Infrared Faucet Control
- Safety Splash Shield
- LED Light Strip
- Magnifier
- Task Light
- Keyboard & Mouse Tray
- Keyboard Arm

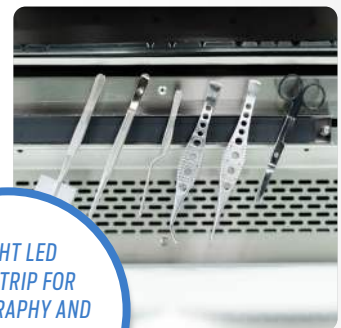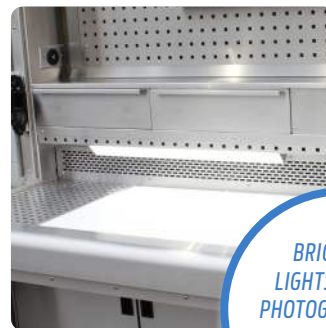

BRIGHT LED LIGHTSTRIP FOR PHOTOGRAPHY AND AN ILLUMINATED WORKSURFACE.

## STANDARD FEATURES

Every Maestro comes with the following included:

- Dimensions: 60" L X 32" W
- Stainless Steel Construction
- Right-Hand Sink
- Backdraft Ventilation
- In-House Ventilation (Dual 8" Duct Stubs)
- Touchscreen Control Panel & Alert System
- On-Demand Front Air System (FAS)
- Flow Detection with Back-Draft Alert Monitoring System
- Seismic Anchoring Points
- Integrated Hood-Mounted Camera Bracket
- Faucet Mixing-Flow Control
- Ensemble Pegboard™ System
- Side Column Accessory Mounts
- 2 GFCI-Protected and 2 USB Receptacles
- 1 Perforated Grid Plate
- Sliding Work Surface White Cutting Board
- 2 Wire Ports for Cord Management
- 1/2 hp Disposal
- Magnetic Instrument Bar
- 1 Large and 1 Small Organizer Bin

# Why choose the Mopec Maestro®?



## Future Flexibility

Should laboratory space or staffing needs change, Maestro benches can be altered within your facility. Mopec technicians have the ability to add pegboard wings or easily replace parts, among other tasks.

## Quality Construction

All Mopec equipment is designed, engineered, and manufactured under the guidelines of ISO 9001:2015 with Design. Every unit is issued with a standard one year warranty.

## Knowledgeable Representatives

Mopec understands that every facility is unique. Our customer service and sales representative teams are educated in common facility challenges and can assist in customizing your equipment.

## Complete Customization

Available in seven different sizes, with over 30 modification options, and a growing inventory of work station accessories, the Mopec Maestro is built to be customizable to your laboratory needs and individual preferences.

## Mopec Maestro® Grossing Stations are ideal for:

Hospital Pathology  
Private Pathology  
Animal Research

Universities  
Histology  
Government Research

Made in the United States | Worldwide Shipping | Over 25 Years of Experience

DESIGN YOUR OWN MOPEC MAESTRO ONLINE AT  
**MOPEC.COM/MAESTRO**



www.mopec.com  
800.362.8491  
800 Tech Row  
Madison Heights, MI 48071