




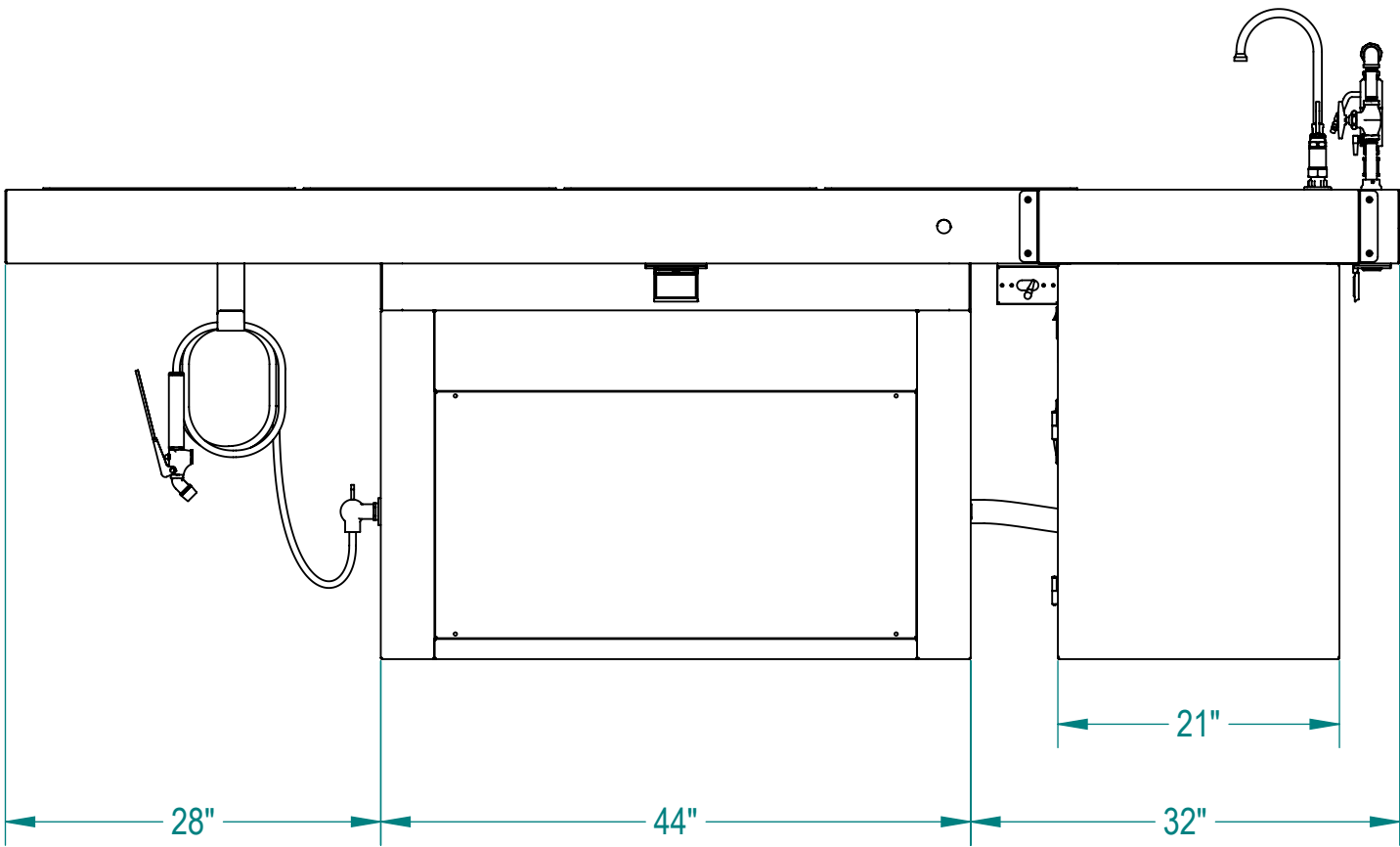

NOTES: UNLESS OTHERWISE SPECIFIED


1. APPLICABLE STANDARDS AND SPECIFICATIONS;
ASME Y14.5-2009, DIMENSIONING AND TOLERANCING
ASME Y14.5-2007, ABBREVIATIONS

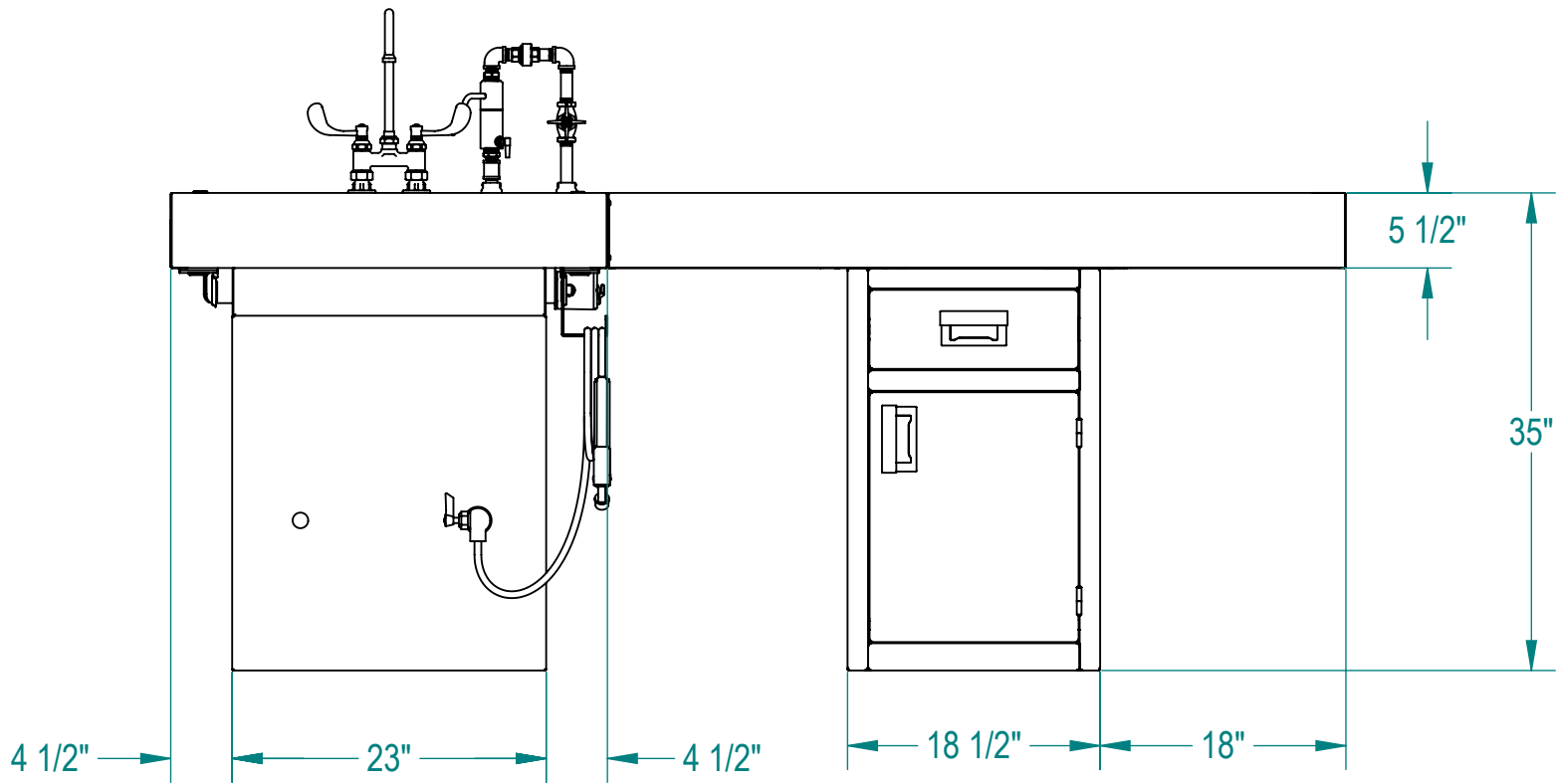

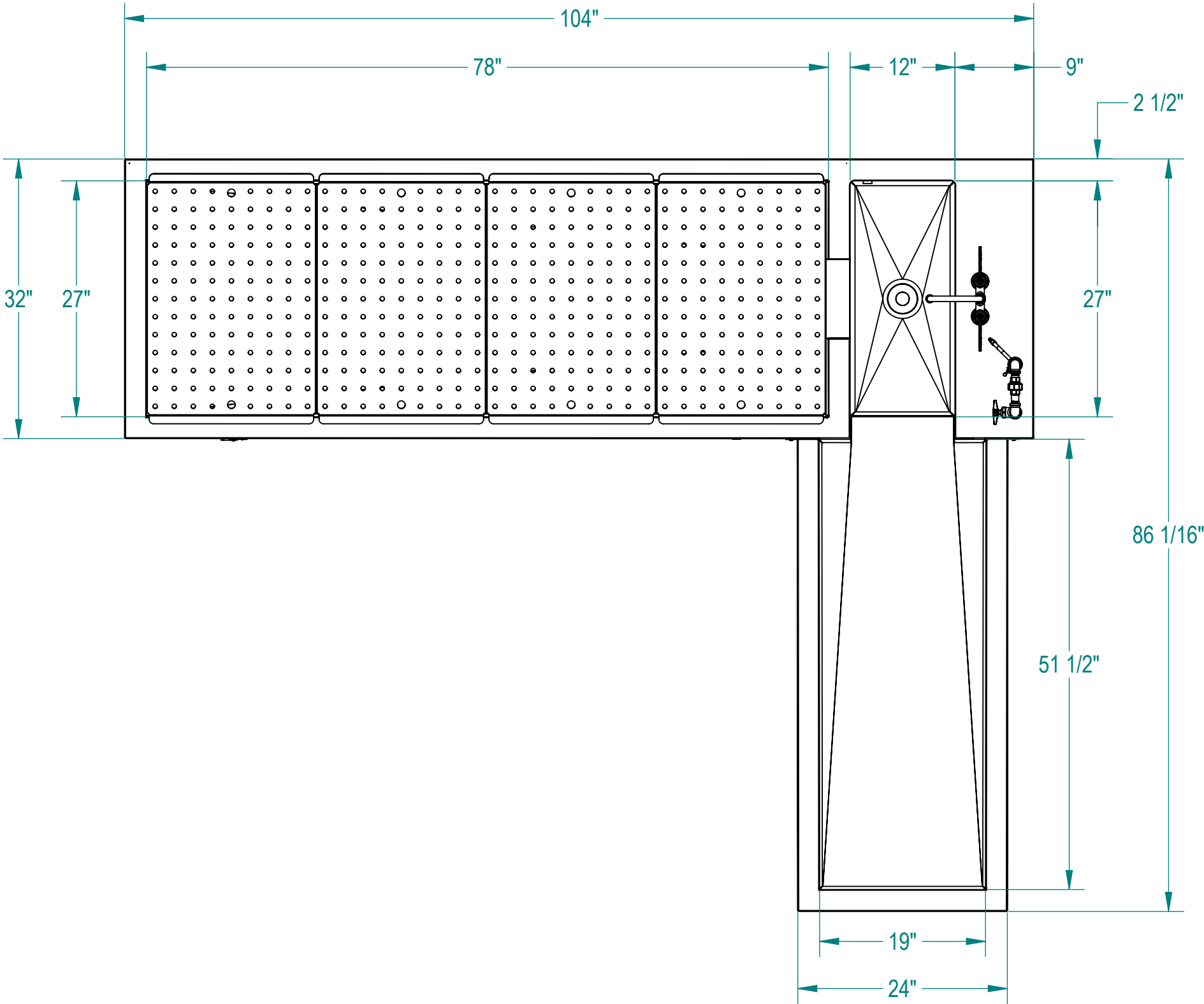

CUSTOMER		APPROVAL		PROPRIETARY THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWIT IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN AUTHORIZATION FROM MOPEC GROUP.						411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.mortechmfg.com	
PROJECT											
				DO NOT SCALE DRAWING		APPROVALS		DATE		TITLE	
				MATERIAL		DRAWN BY		11/05/24		MODEL 1036-19R (CE600)	
				304 STAINLESS STEEL		ELISEO S.				PEDESTAL AUTOPSY TABLE, L-SHAPED, RIGHT	
				FINISH		CHECKED BY				SIZE	
				#4		APPROVED BY				DWG NO.	
										1036-19R	
										REV.	
										A	
										SCALE	
										1:14	
										NEXT ASSY.	
										SHEET	
										1 OF 4	

NOTES: UNLESS OTHERWISE SPECIFIED

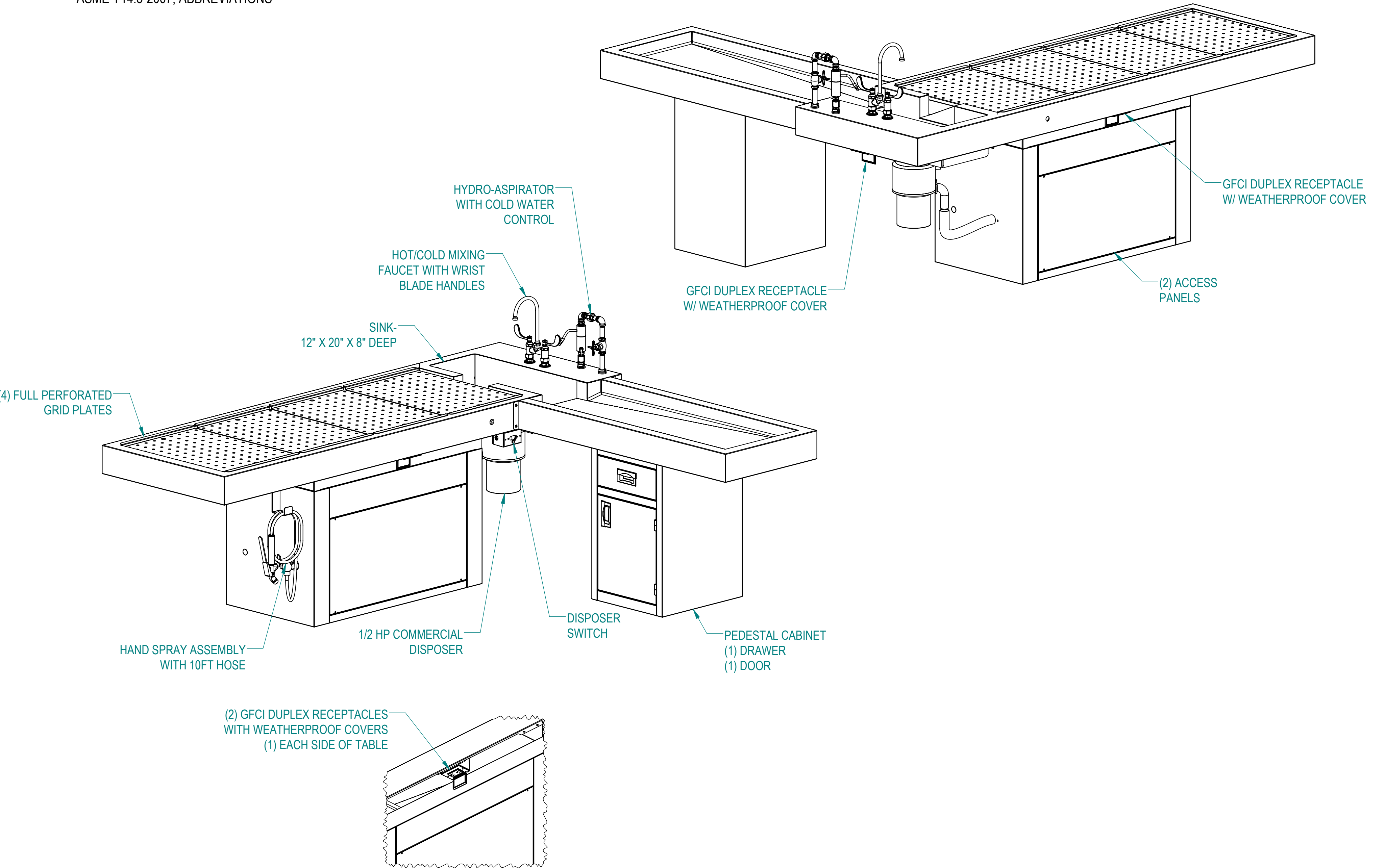
1. APPLICABLE STANDARDS AND SPECIFICATIONS;
ASME Y14.5-2009, DIMENSIONING AND TOLERANCING
ASME Y14.5-2007, ABBREVIATIONS



CUSTOMER		APPROVAL		PROPRIETARY THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWIT IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN AUTHORIZATION FROM MOPEC GROUP.						411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.mopectechmfg.com	
PROJECT						TITLE					
				DO NOT SCALE DRAWING		APPROVALS		DATE		MODEL 1036-19R (CE600) PEDESTAL AUTOPSY TABLE, L-SHAPED, RIGHT	
				MATERIAL		DRAWN BY		11/05/24			
				304 STAINLESS STEEL		ELISEO S.					
				FINISH		CHECKED BY				REV.	
				#4		APPROVED BY				A	
								SCALE		1:12	
								NEXT ASSY.			
								SHEET		2 OF 4	

UTILITY REQUIREMENTS:

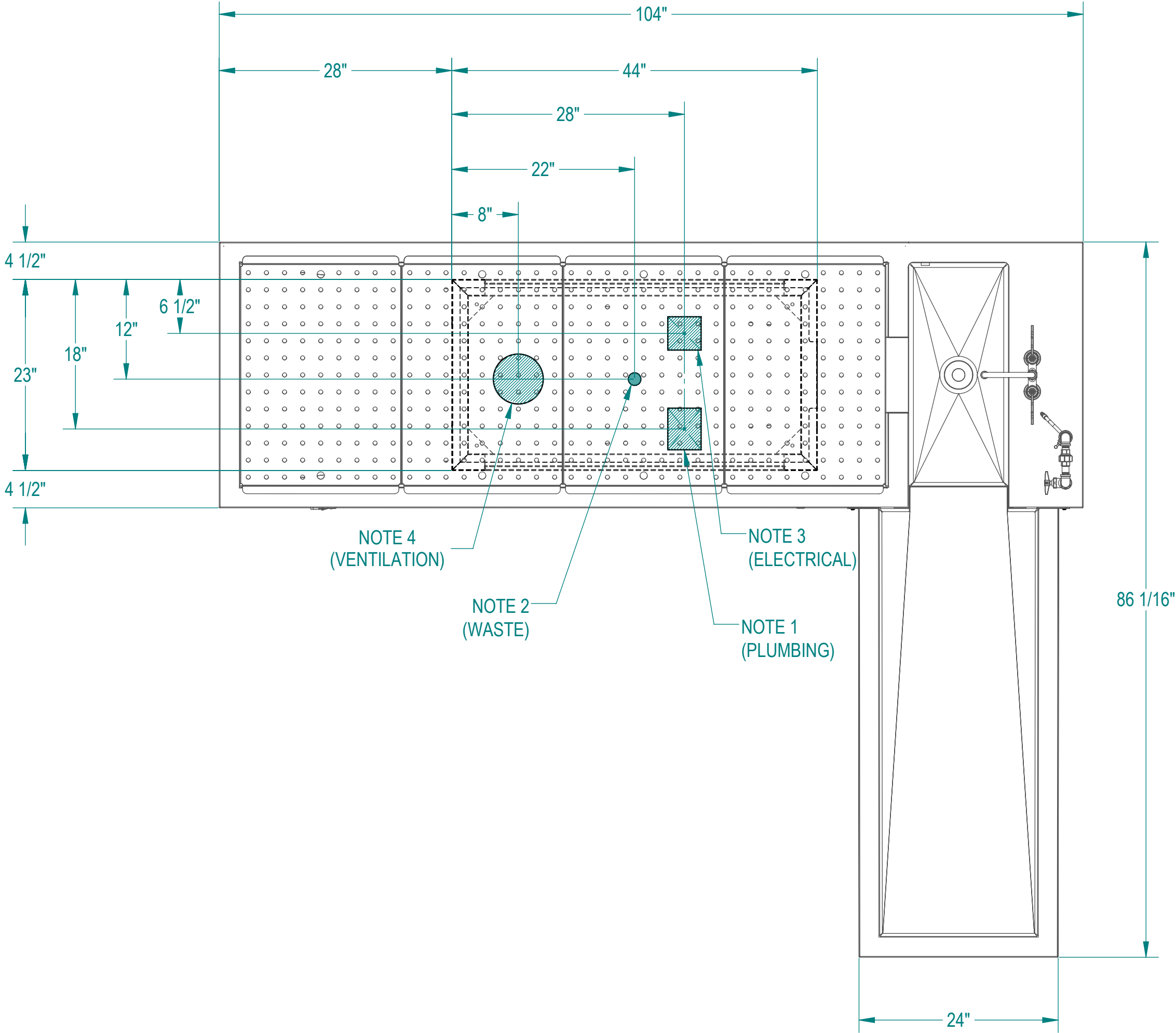

1. A 4" X 5" AREA FOR INLET PLUMBING TO INCLUDE THE FOLLOWING:
(1) COLD WATER SUPPLY, 3/4" OD COPPER STUBBED OUT 6"
(1) HOT WATER SUPPLY, 1/2" OD COPPER STUBBED OUT 6"
SUPPLY SHUT-OFF VALVES AT EACH WATER SUPPLY
FLUSH LINES CLEAR BEFORE CONNECTING TO STATION
2. (1) STANDARD 1 1/2" WASTE LINE STUBBED OUT 6"
SUPPLY P-TRAP FOR DRAIN
3. A 4" X 4" AREA FOR ELECTRICAL INPUT TO INCLUDE THE FOLLOWING:
(1) 115V/60Hz/1PH/15A CIRCUIT FOR STATION
(1)115V/60Hz/1PH/15A CIRCUIT FOR DISPOSER
4. A 6" DIA DUCT EXHAUST CONNECTION:
RECOMMENDED EXHAUST: 450-500 CFM AT 0.50 STATIC PRESSURE
MEASURED AT THE 6" EXHAUST, PROVIDE WATER CAP OVER DUCT INLET.
EXISTING CONDITIONS MAY WARRANT AN INCREASE OR DECREASE,
CONTACT YOUR ENGINEERING DEPT.
5. PROVIDE 60-65 PSI WATER PRESSURE AND 4 GALLONS PER MINUTE OF
FLOW FOR PROPER OPERATION OF HYDRO ASPIRATOR
6. ELECTRICIAN TO PROVIDE LIQUID TIGHT JUNCTION BOX AND CONDUIT PER
LOCAL CODES.



CUSTOMER	APPROVAL		<div>PROPRIETARY</div> <div>THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWITH IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN AUTHORIZATION FROM MOPEC GROUP.</div>			<div><div><div>MOPEC</div><div>GROUP</div></div><div>411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.mortechmfg.com</div></div>			
PROJECT						<div>TITLE</div> <div>MODEL 1036-19R (CE600)</div> <div>PEDESTAL AUTOPSY TABLE, L-SHAPED, RIGHT</div>			
			DO NOT SCALE DRAWING		APPROVALS	DATE	SIZE DWG NO. REV.		
			MATERIAL	DRAWN BY	11/05/24		C 1036-19R A		
			304 STAINLESS STEEL	ELISEO S.					
			FINISH	CHECKED BY					
			#4		APPROVED BY		SCALE 1:12 NEXT ASSY. SHEET 3 OF 4		

1. APPLICABLE STANDARDS AND SPECIFICATIONS;
ASME Y14.5-2009, DIMENSIONING AND TOLERANCING
ASME Y14.5-2007, ABBREVIATIONS



CUSTOMER	APPROVAL		<p align="center"><u>PROPRIETARY</u></p> <p align="center">THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWIT IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN AUTHORIZATION FROM MOPEC GROUP.</p>			  <div> 411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.mopecsmfg.com </div>	
PROJECT						TITLE	
						<p align="center">MODEL 1036-19R (CE600)</p> <p align="center">PEDESTAL AUTOPSY TABLE, L-SHAPED, RIGHT</p>	
						SIZE	REV.
						C	A
	SIGNATURE	DATE	DO NOT SCALE DRAWING	APPROVALS	DATE	SCALE	
			MATERIAL	DRAWN BY		DWG NO.	
			304 STAINLESS STEEL	ELISEO S.	11/05/24	1036-19R	
			FINISH	CHECKED BY		NEXT ASSY.	
			#4	APPROVED BY		SCALE	
						1:12	
						SHEET	
						4	
						OF	
						4	