


Maestro Hand Spray and End Rinse



MB1036R

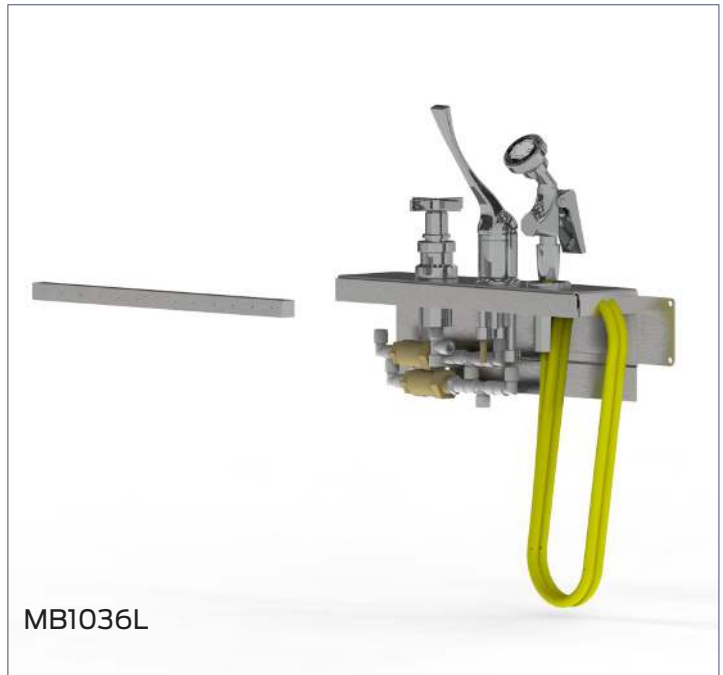

MB1036L

SPECIFICATIONS

Construction

- Designed for customers to rinse off the entire workstation of their grossing station
- Includes (1) sprayer, (1) 60" hose, (1) end rinse bar, flow control fixture, plumbing and fittings
- The flow for the hand spray and for the end rinse are controlled by separate valves

Dimensions

- 60" hose
- Rinse bar: 16.3" L x .75" W x .75" H

Materials

- Hand sprayer with soft spray: chrome-plated brass
- Rinse Bar and Sheetmetal: 304 stainless steel construction with a polished No. 4 finish
- Hardware: 304 stainless steel
- Plumbing lines and fittings: brass, pex, and stainless steel

Performance

- Mixing valve/flow control fixture with vacuum breaker
- Backflow preventers on both hot and cold-water lines

Options & Accessories

- MB1036L: designed for left-hand Maestro grossing stations
- MB1036R: designed for right-hand Maestro grossing stations

Compatibility

- Compatible with all Maestro grossing stations that do not currently have this option

NOTE:

- For stations adding this option where it had not previously existed, it would include a new endcap to the worksurface with a hole for the hand spray to retract inside the station
- The Hand Spray Rinse **only** now comes standard in all Maestro Encore grossing stations

Variations

Model No.	Description	Compatibility
MB1036L	Hand Spray & End Rinse, Left-Hand Units	Left-Hand Maestro Grossing Stations
MB1036R	Hand Spray & End Rinse, Right-Hand Units	Right-Hand Maestro Grossing Stations

MOPEC
GROUP

mopec
ELEVATING PATHOLOGY

MORTECH
MANUFACTURING

www.mopecgroup.com

800.362.8491

800 Tech Row, Madison Heights, MI