



# Mopec Maestro Encore Ventilation Options

## In-House Ventilation FACILITY-PROVIDED, STANDARD OFFERING

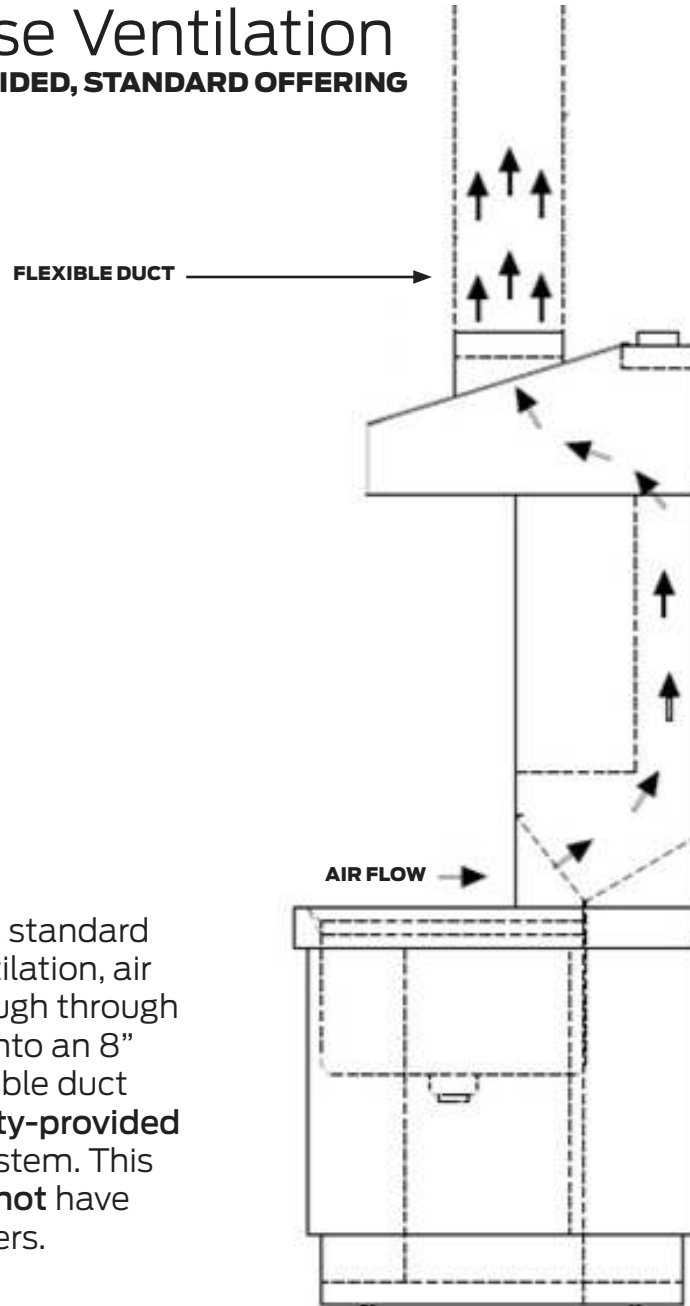

With Mopec's standard in-house ventilation, air is drawn through through the unit and into an 8" diameter flexible duct into the **facility-provided** ventilation system. This system does **not** have internal blowers.

## Recirculating Option PURCHASED AT ADDITIONAL COST

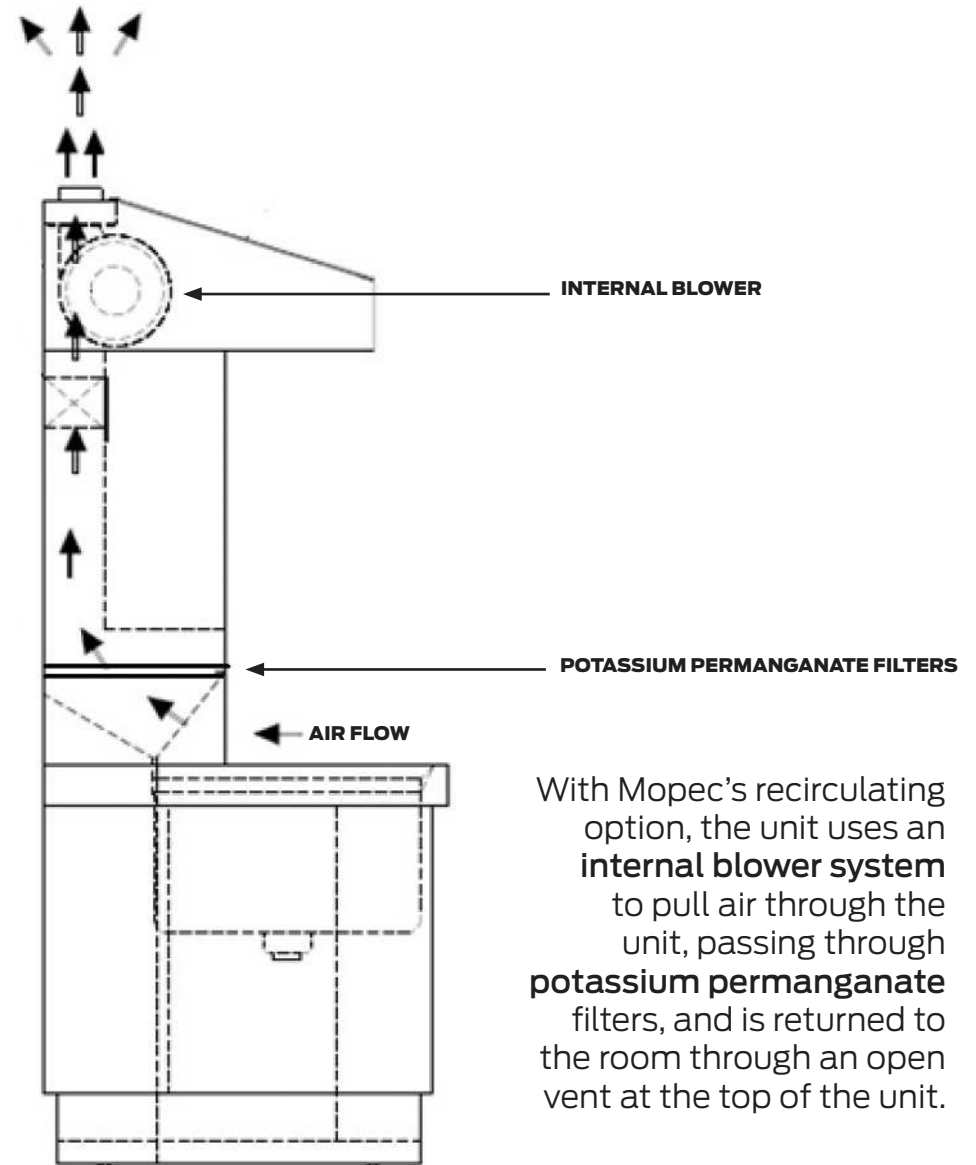

With Mopec's recirculating option, the unit uses an **internal blower system** to pull air through the unit, passing through **potassium permanganate filters**, and is returned to the room through an open vent at the top of the unit.