
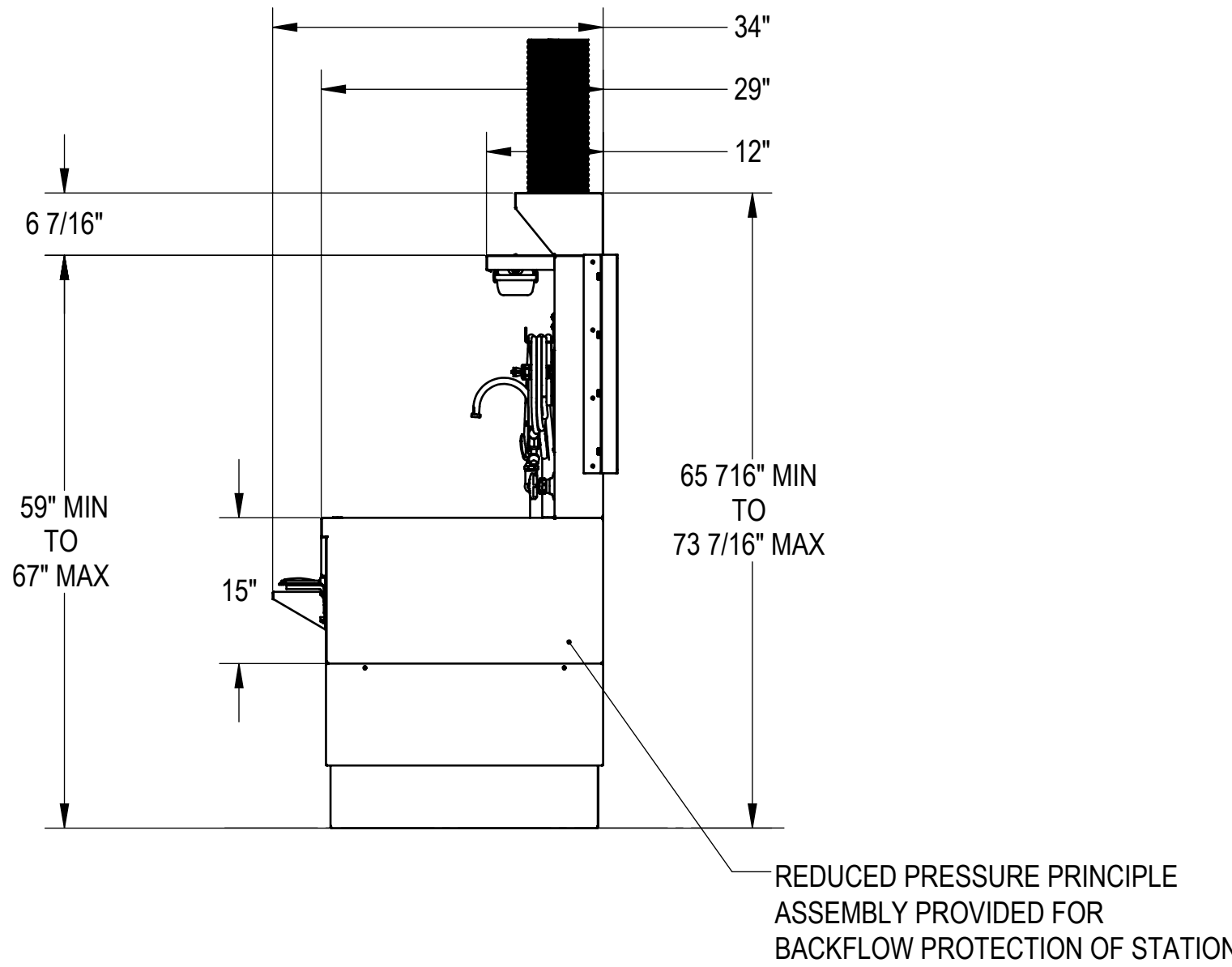
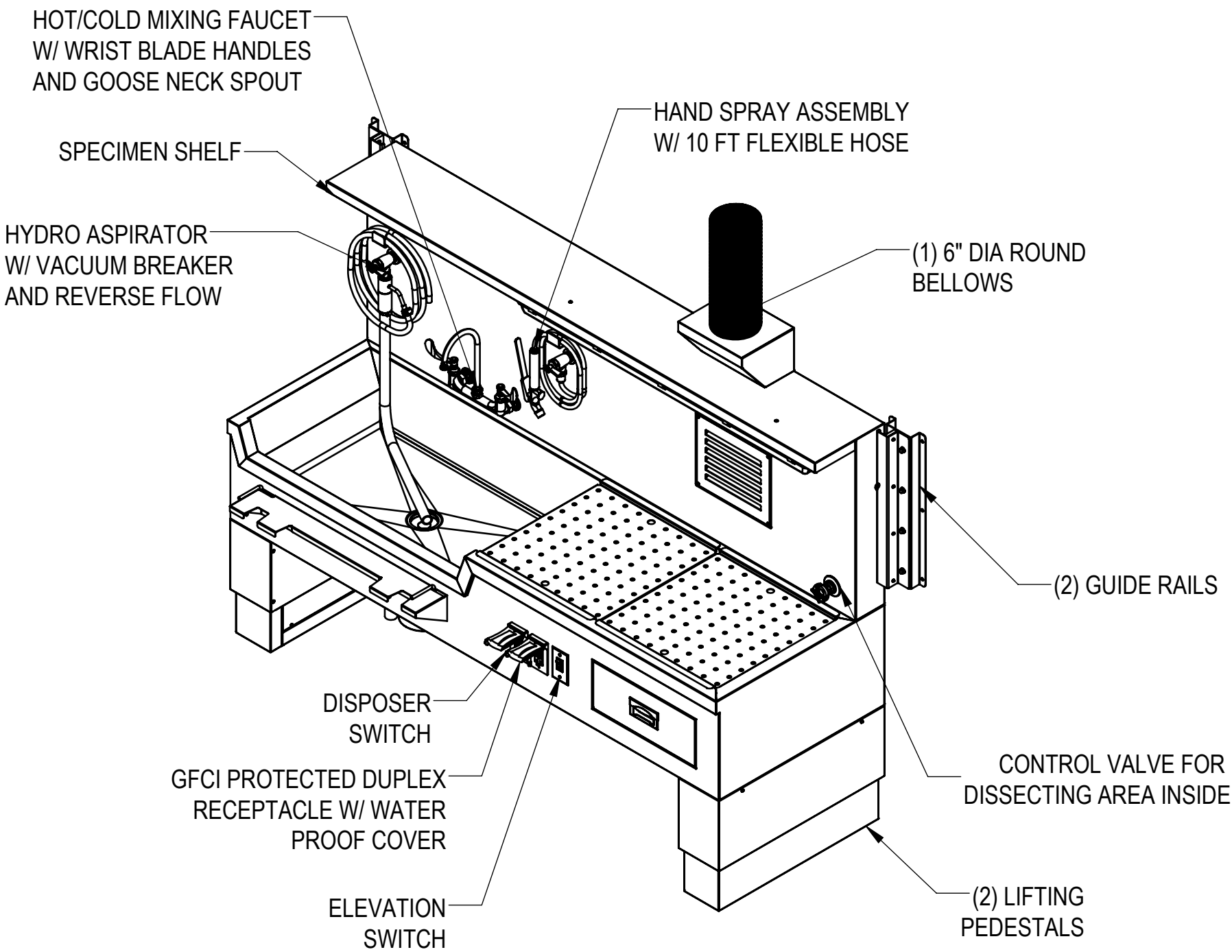


NOTES: UNLESS OTHERWISE SPECIFIED



1. APPLICABLE STANDARDS AND SPECIFICATIONS;
ASMEY14.5-2009, DIMENSIONING AND TOLERANCING
ASME Y14.5-2007, ABBREVIATIONS

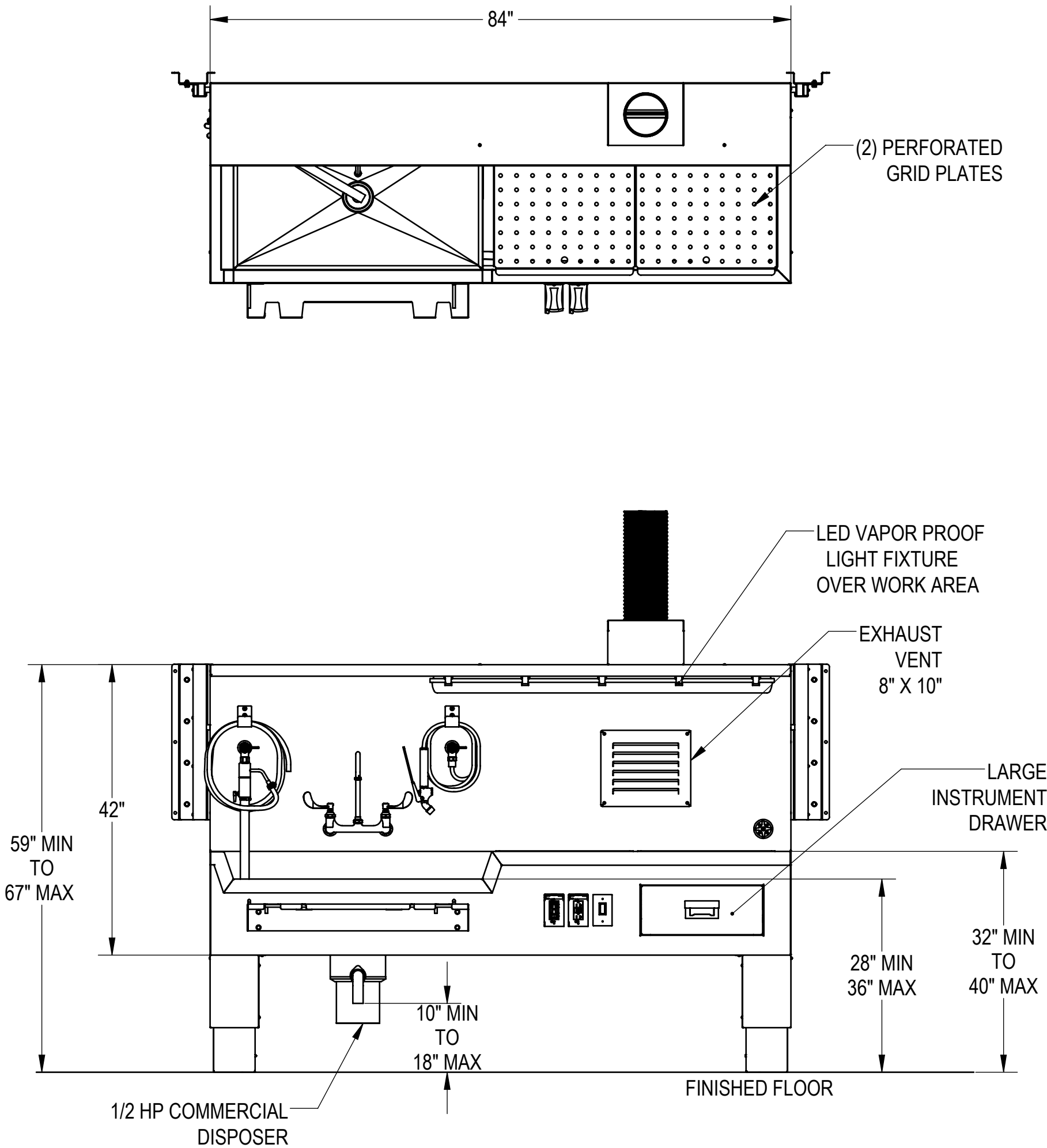

CUSTOMER		APPROVAL		<div>PROPRIETARY</div> <div>THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWITH IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN AUTHORIZATION FROM MORTECH MFG.</div>			<div><div><div>MOPEC GROUP</div><div><div>MORTECH</div><div><small>MORTECH GROUP</small></div></div></div><div>ISO 9001:2015 CERTIFIED COMPANY</div><div>411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.mortechmfg.com</div></div>					
PROJECT							TITLE					
							MODEL 1036-7A FLOOR MOUNT AUTOPSY STATION, LEFT SINK, ELEVATING					
		SIGNATURE		DATE		SIZE			DWG NO.		REV.	
						C			1036-7A-AP		X01	
						SCALE			1:16		NEXT ASSY.	
									SHEET		1 OF 4	

UTILITY REQUIREMENT NOTES



1. FLUSH LINES CLEAR BEFORE CONNECTING TO STATION.
2. PLUMBER TO PROVIDE SHUT OFF VALVES FOR WATER CONNECTIONS AND P-TRAP FOR DRAIN.
3. PROVIDE 60-65 PSI WATER PRESSURE AND 4 GALLONS PER MINUTE OF FLOW FOR PROPER OPERATION OF HYDRO ASPIRATOR.
4. ELECTRICIAN TO PROVIDE LIQUID TIGHT JUNCTION BOX AND CONDUIT CONNECTIONS PER LOCAL CODES.
5. RECOMMENDED EXHAUST: 300-350 CFM PER EXHAUST VENT WITH APPROXIMATE PRESSURE DROP OF 0.25 W.G.

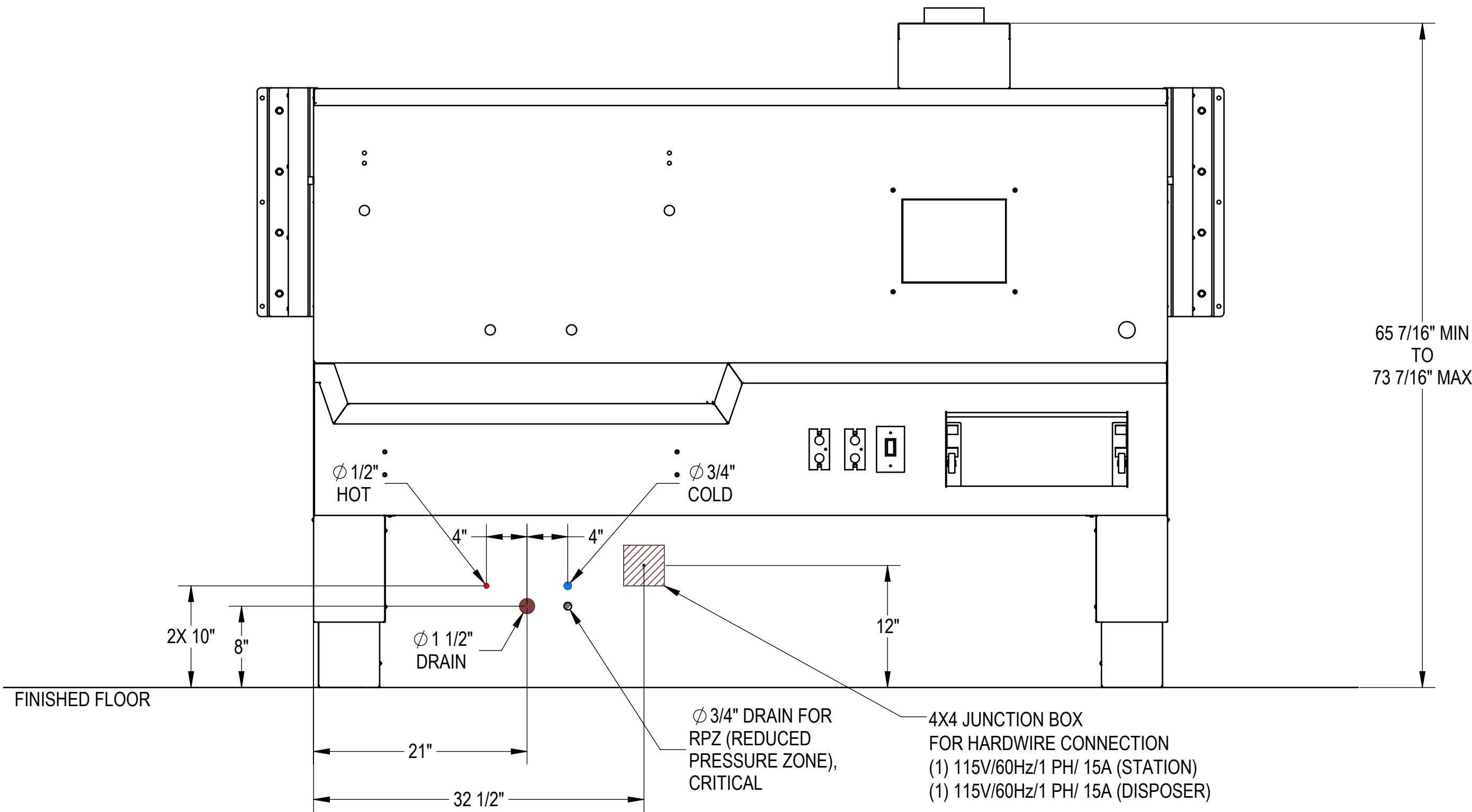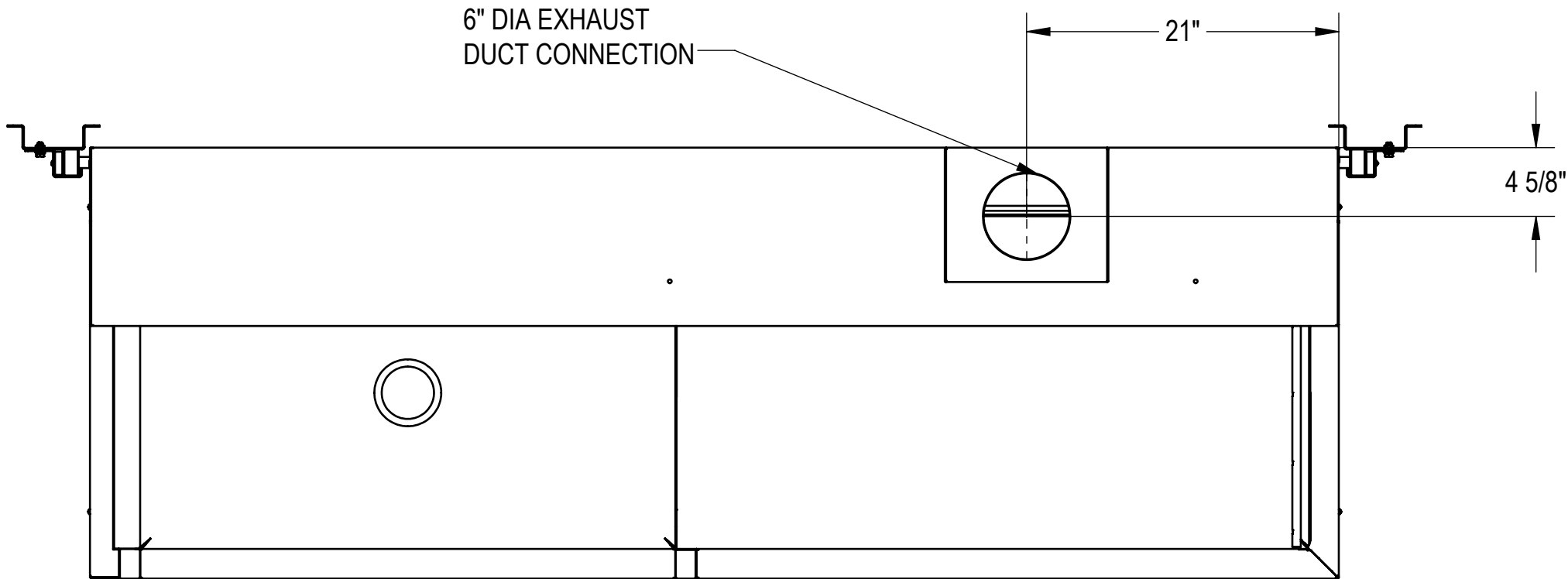

CUSTOMER	APPROVAL		<div>PROPRIETARY</div> <div>THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWITH IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN AUTHORIZATION FROM MORTECH MFG.</div>			<div><div><div>MOPEC GROUP</div><div>MopecMORTECH</div></div><div>ISO 9001:2015 CERTIFIED COMPANY</div><div>411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.mortechmfg.com</div></div>					
PROJECT						<div>TITLE</div> <div>MODEL 1036-7A</div> <div>FLOOR MOUNT AUTOPSY STATION, LEFT SINK, ELEVATING, ROUGH-IN</div>					
SIGNATURE			DATE			SIZE		DWG NO.		REV.	
						C		1036-7A-AP		X01	
						SCALE 1:10		NEXT ASSY.		SHEET 2 OF 4	
			DO NOT SCALE DRAWING			APPROVALS		DATE			
			MATERIAL			DRAWN BY					
			304 STAINLESS STEEL			LUIS F.		06/06/19			
			FINISH			CHECKED BY					
			#4 SATIN			APPROVED BY					

NOTES: UNLESS OTHERWISE SPECIFIED


1. INSTALLATION IS BASED ON ROUGH-IN MOUNTING
DETAIL PROVIDED.
2. GUIDE RAILS PREVENT TOP OF SINK FROM
PULLING AWAY FROM WALL.

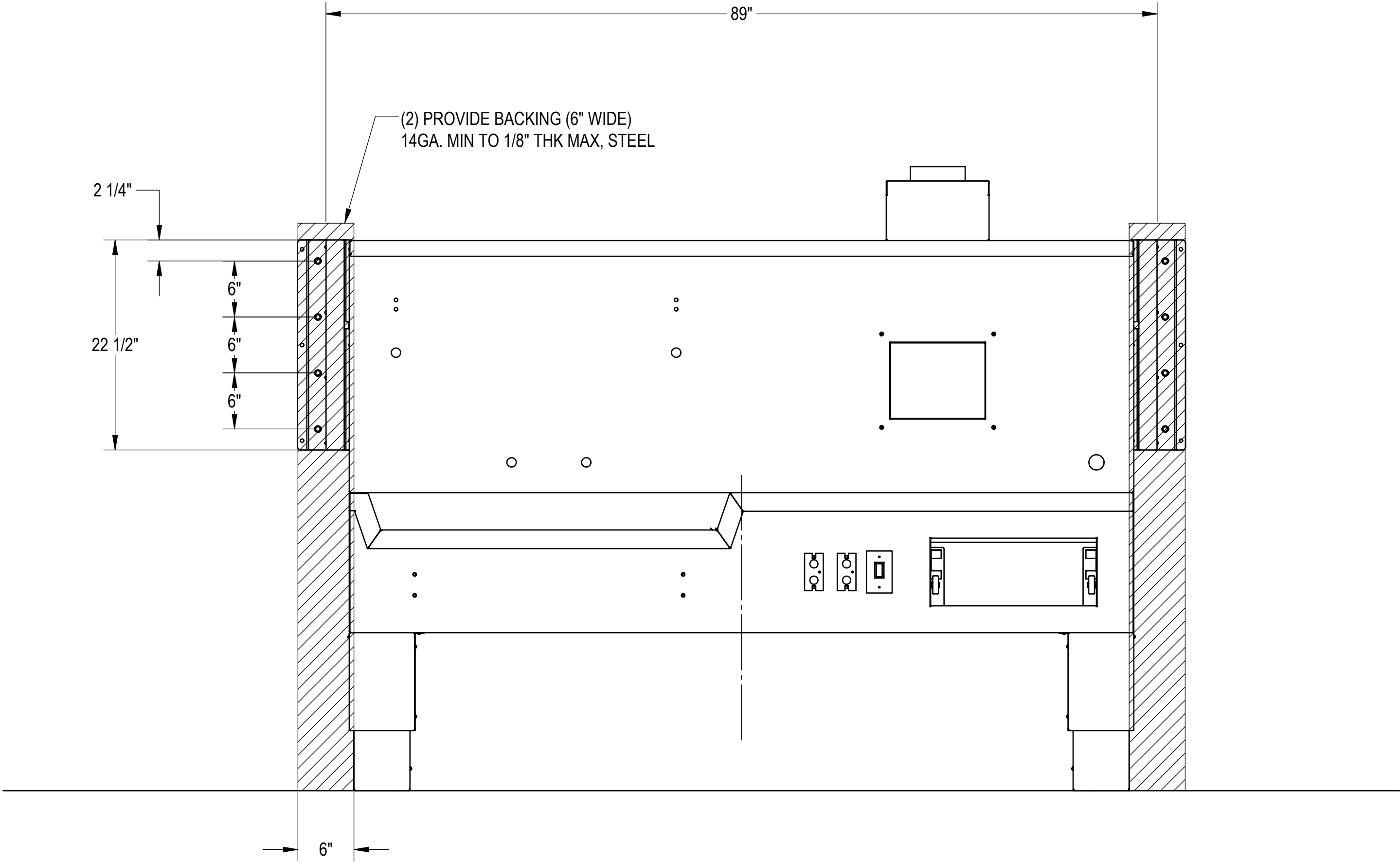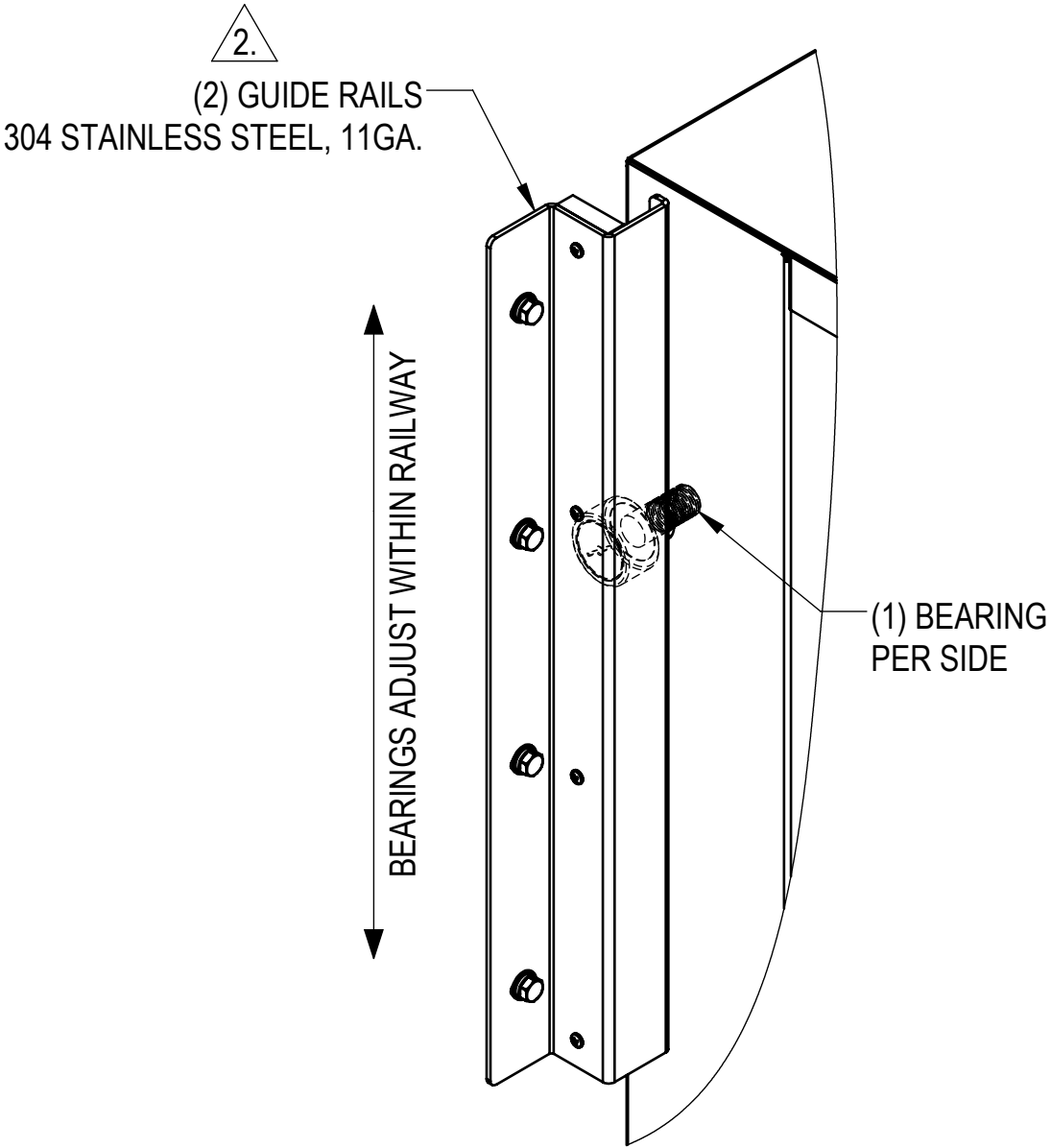

CUSTOMER	APPROVAL		<div>PROPRIETARY</div> <div>THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWITH IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN AUTHORIZATION FROM MORTECH MFG.</div>			<div><div><div>MOPEC GROUP</div><div>MORTECH</div></div><div>ISO 9001:2015 CERTIFIED COMPANY</div><div>411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.mortechmfg.com</div></div>					
PROJECT											
			DO NOT SCALE DRAWING			APPROVALS		DATE			
			MATERIAL		DRAWN BY		LUIS F.		06/06/19		
			304 STAINLESS STEEL		CHECKED BY						
			FINISH		APPROVED BY						
			#4 SATIN								
						SCALE		1:10		NEXT ASSY.	
						SHEET		3		OF 4	
						TITLE		MODEL 1036-7A FLOOR MOUNT AUTOPSY STATION, LEFT SINK, ELEVATING, WALL RAIL			
						SIZE		DWG NO.		REV.	
						C		1036-7A-AP		X01	



NOTES: UNLESS OTHERWISE SPECIFIED

1. APPLICABLE STANDARDS AND SPECIFICATIONS;
ASMEY14.5-2009, DIMENSIONING AND TOLERANCING
ASME Y14.5-2007, ABBREVIATIONS

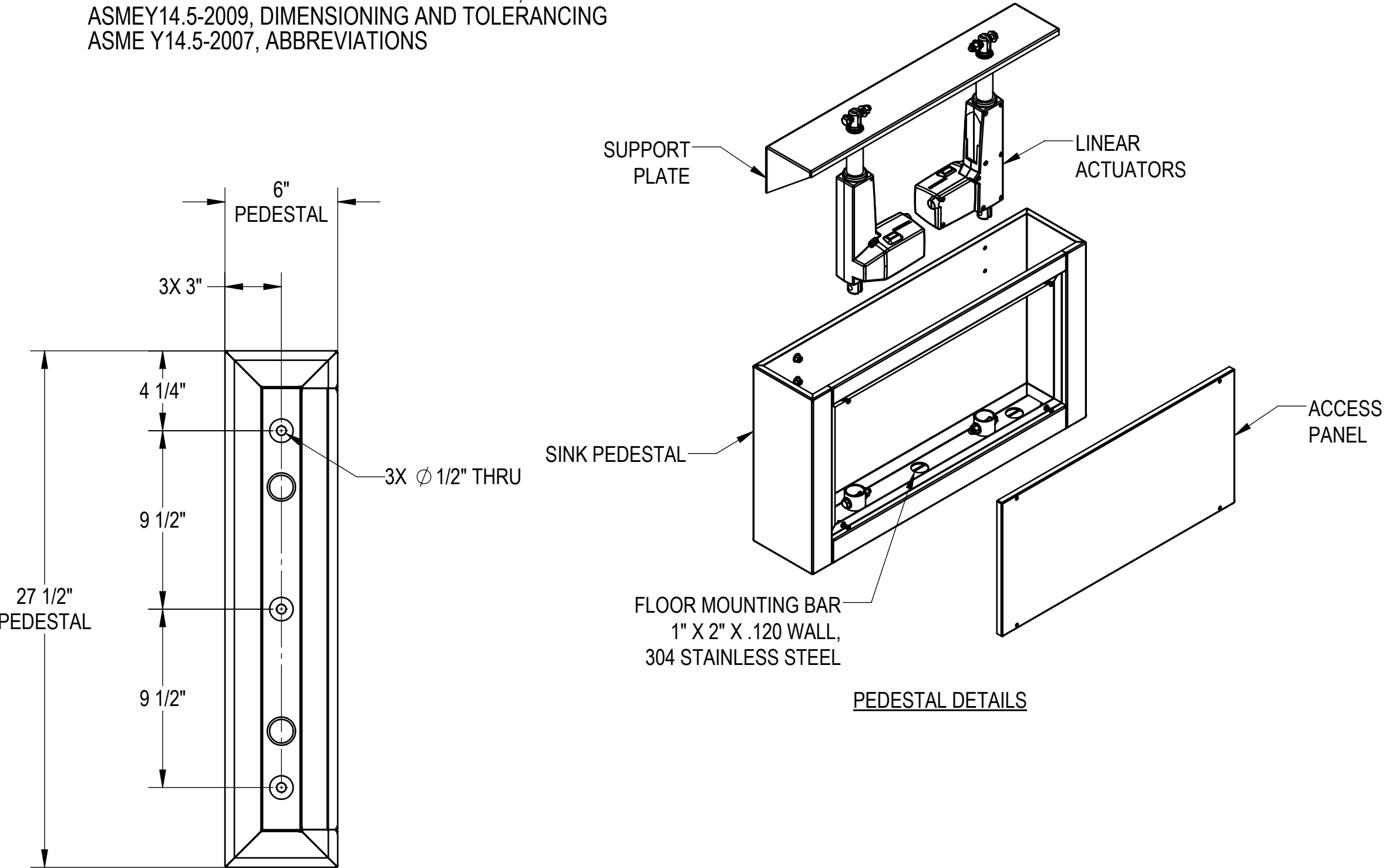

PEDESTAL DETAILS

PEDESTAL TOP VIEW
(2) PER UNIT

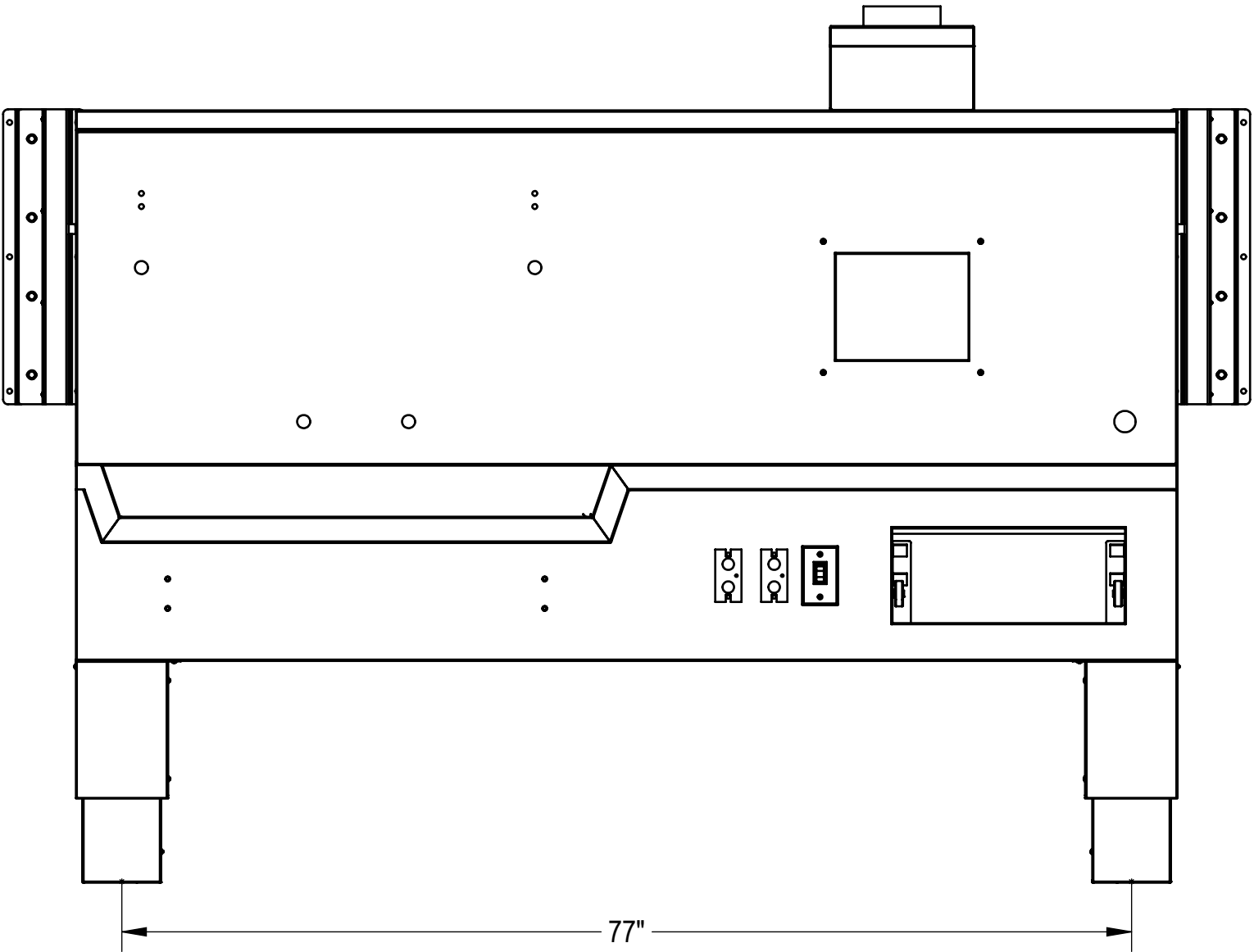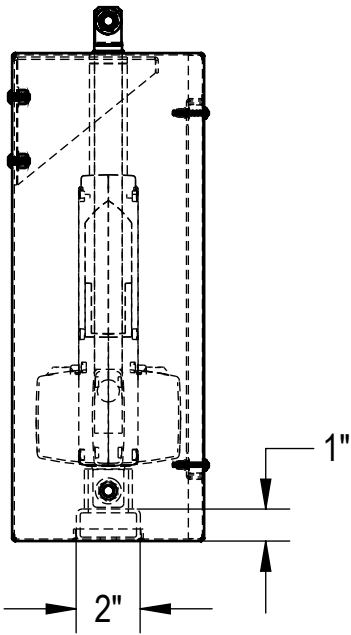

CUSTOMER	APPROVAL		PROPRIETARY			<div>MOPEC GROUP</div> <div>Mopec Inc. MORTECH</div> <div>ISO 9001:2015 CERTIFIED COMPANY</div> <div>411 North Aerojet Avenue Azusa, CA 91702 tel: (800) 410-0100 (626) 334-1471 fax: (626) 334-1704 www.moritechmfg.com</div>			
PROJECT			THIS DOCUMENT AND THE DATE DISCLOSED HEREIN AND HEREWITH IS NOT TO BE REPRODUCED, USED, OR OTHERWISE DISCLOSED IN WHOLE OR IN PART TO ANYONE WITHOUT WRITTEN AUTHORIZATION FROM MORTECH MFG.			TITLE			
						MODEL 1036-7A			
						FLOOR MOUNT AUTOPSY STATION, LEFT SINK, ELEVATING FLOOR MOUNTING			
	SIGNATURE		DATE	DO NOT SCALE DRAWING	APPROVALS	DATE	SIZE	DWG NO.	REV.
			MATERIAL	DRAWN BY		C			
			304 STAINLESS STEEL	LUIS F.	06/06/19				
			FINISH	CHECKED BY					
			#4 SATIN	APPROVED BY		SCALE	1:10	NEXT ASSY.	
								SHEET	
								4 OF 4	