


Maestro Dry Contact for Building Management System

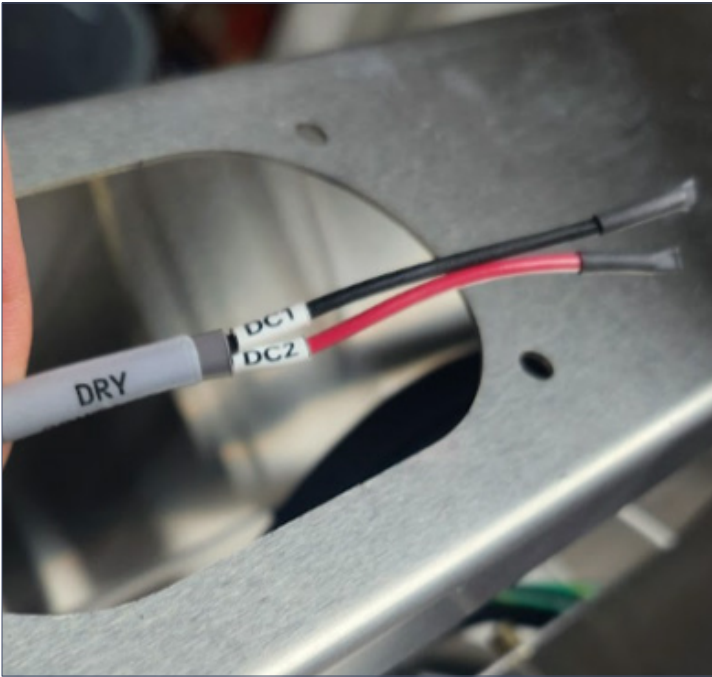

SPECIFICATIONS

Construction

- Electrical wires to connect a Maestro or Maestro Encore grossing station to the facility Building Management System
- Signal wires run from facility to the grossing station and connect to the alarm relay signal output on electrical cabinet DIN rail
- Includes one set of dry contacts

Performance

- This feature connects to the facility's Building Management System to inform the appropriate parties that the grossing station has an alert condition
- BMS output triggering is setup in the Notifications Setting Page in the touchscreen panel

Compatibility

- Compatible with all Maestro and Maestro Encore grossing stations
- Maestro Encore units have this option integrated into the wiring harness and labeled "DRY CONTACTS DC1/DC2)"

NOTE:

- We recommend ordering at the time of purchase to be properly installed during manufacturing of the grossing station.

MOPEC
GROUP

mopec
ELEVATING PATHOLOGY

MORTECH
MANUFACTURING

www.mopecgroup.com

800.362.8491

800 Tech Row, Madison Heights, MI