



Autopsy Station, Center Sink

Mopec's Autopsy Stations are a durable, stainless steel workstation designed for professional pathology environments. It offers a spacious, easy-to-clean work surface with efficient drainage, integrated utilities, and features that support smooth workflow and sanitation. With built-in lighting, storage, waste handling, and compatibility with Mopec carts, this line is reliable for everyday use in demanding settings.

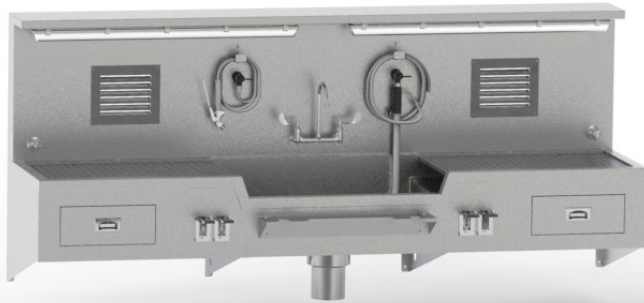

1036-10

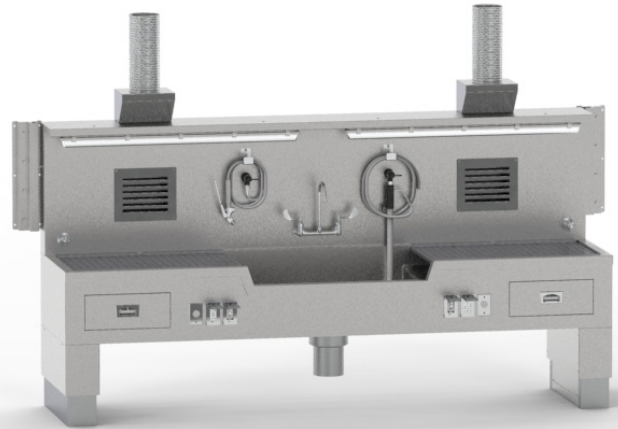

1036-10A

SPECIFICATIONS

Construction

- 1036-10: Wall Mount
- 1036-10A: Floor Mount, 8" Elevation
- Grid plate work surfaces with 1/2" recessed holes to allow for drainage
- 10' Spray hose assembly that can reach the length of the sink
- (2) Large instrument drawers
- (2) End rinses that drain into a common sink
- Cart latch works with most Mopec Carts except those with hydraulic bases
- Wall mounted passive ventilation through Facility-provided HVAC
- Cart-latch mechanism with operation handle
- (2) GFCI duplex receptacle with waterproof cover
- Heavy-duty commercial waste disposer, 1/2 hp and 2 switches
- Light fixtures over the work areas
- Deluxe aspirator with reverse flow to relieve suction clogging

Utilities / Power

- 2-duplex GFCI 120V / 60Hz outlets, Requires 2- 15-amp circuits, one for station, one for disposal
- 1/2" Hot water, 3/4" Cold Water connection dropping to 1/2" after RPZ valve, and 1-1/2" Drain for water hookup
- Provide 60-65 psi of water pressure at 4 Gallons per minute for proper operation of hydro aspirator
- Uses facility ventilation Capable of 700CFM at 0.25 W.G
- Recommended Supply 350 CFM X 2 Grilles, 700 CFM total

Materials

- 14-18 gauge 304 Stainless steel with No. 4 finish
- Pex piping for plumbing purposes
- Chrome Plated Brass fixtures

Dimensions

- 29" L X 114" W (123" including travel brackets) X 50" H
- Sink Dimensions: 16" L X 30" W X 6" D

Variations

Model No.	Description	Dimensions
1036-10	Autopsy Station, Center Sink, Wall Mount	29"L x 114"W x 50"H
1036-10A	Autopsy Station, Center Sink, Floor Mount	29"L x 114"W x 50"H

MOPEC
GROUP

mopec
ELEVATING PATHOLOGY

MORTECH
MANUFACTURING

www.mopecgroup.com

800.362.8491

800 Tech Row, Madison Heights, MI