


Autopsy Station, Right Sink

Mopec's Autopsy Stations are a durable, stainless steel workstation designed for professional pathology environments. It offers a spacious, easy-to-clean work surface with efficient drainage, integrated utilities, and features that support smooth workflow and sanitation. With built-in lighting, storage, waste handling, and compatibility with Mopec carts, this line is reliable for everyday use in demanding settings.

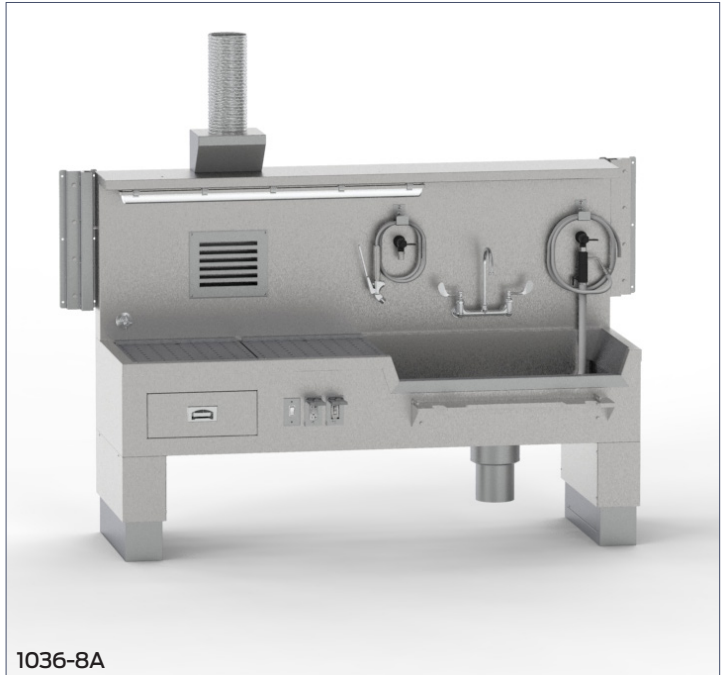

SPECIFICATIONS

Construction

- TÜV Certified
- 1036-8: Wall Mount
- 1036-8A: Floor Mount, 8" elevation with electromechanical linear actuator
- Hot / cold water fixture with gooseneck faucet and wrist blade handles
- Exhaust duct with 8" x 10" grill
- Dissecting area rinse assembly with cold water valve
- Heavy-duty commercial waste disposer, 1/2 HP
- Deluxe aspirator with reverse flow to relieve suction clogging
- GFCI duplex receptacles with waterproof cover
- Flexible 10 ft. spray hose, with cold water control valve and nozzle
- Full length specimen shelf (12" W)
- Large instrument drawer
- Light fixtures over work area
- Cart-to-sink locking mechanism
- Sink on right side

Utilities / Power

- Water: Provide 60-65 PSI water pressure and 4 GPM of flow
Refer to rough-in drawing for plumbing connection sizes
- Exhaust: 600-650 CFM per grille with approximate pressure drop of 0.25 W.G
- Electrical: Provide (2) circuits in 4" x 4" junction box:
(1) 115V/60Hz / 1Ph / 15A (for station)
(1) 115V/60Hz / 1Ph / 15A (for disposal)

Materials

- 304 Stainless steel with No. 4 finish

Dimensions

- 34" L X 84" W
- 1036-7: 50" H or 61" H (mounted)
- 1036-7A: 59" - 67" H (8" elevation)
- Note: Height excludes exhaust connection
- Sink: 20.75" L X 35" W X 6" D

Variations

Model No.	Description	Dimensions
1036-8	Autopsy Station, Right Sink, Wall Mount	34"L x 84"W x 61"H
1036-8A	Autopsy Station, Right Sink, Floor Mount	34"L x 84"W x 67"H

MOPEC
GROUP

mopec
ELEVATING PATHOLOGY

MORTECH
MANUFACTURING

www.mopecgroup.com

800.362.8491

800 Tech Row, Madison Heights, MI